



| | | |
|--------|---------------------------|--------------|
| Title | | |
| CURTIN | | |
| Size | Document Number | Rev. |
| A | <Doc> | <Rev Code> |
| Date: | Friday, February 01, 2008 | Sheet 1 of 1 |
| | 2 | 1 |

THE INFORMATION CONTAINED HEREIN IS PROPRIETARY TO COMPU-MEDICS AND IS TO BE USED BY THE RECIPIENT SOLELY FOR THE PURPOSE FOR WHICH IT IS SUPPLIED. IT SHALL NOT BE DUPLICATED, IN WHOLE OR IN PART WITHOUT PRIOR WRITTEN PERMISSION FROM AN AUTHORIZED REPRESENTATIVE OF COMPU-MEDICS.


| REVISIONS | | | | |
|-----------|-----------------------------------|----------|------|----------|
| REV | DESCRIPTION | DCN | INIT | DATE |
| A | INITIAL RELEASE | 0401-001 | WC | 01/05/04 |
| 8218A | Correction to allow drawing as PN | 0403-002 | WC | 03/02/04 |

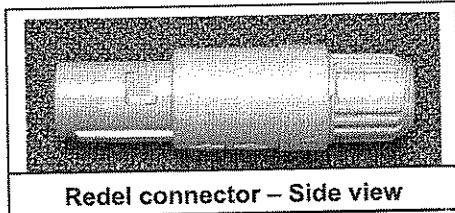

Redel connector – Side view

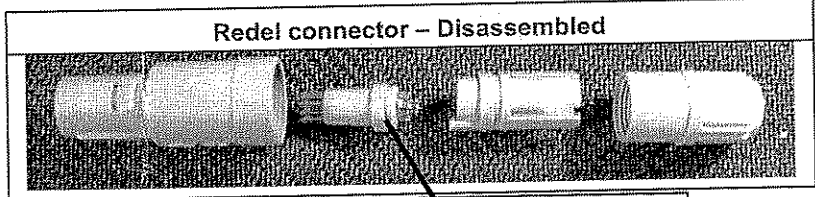

Redel connector – Disassembled

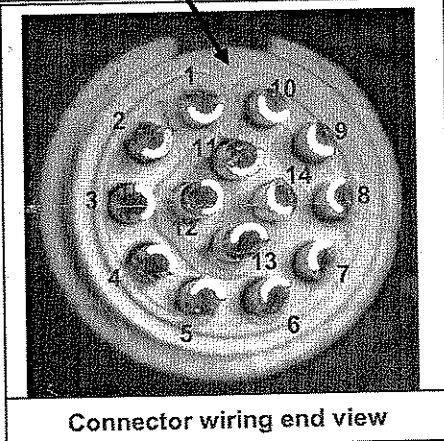

Connector wiring end view

REDEL HIGH LEVEL PINOUT

| PIN | INPUT/ OUTPUT | FUNCTION |
|-----|---------------|----------------|
| 1 | Input | Hi-Level 1 (+) |
| 2 | Output | + DC Voltage |
| 3 | Ground | |
| 4 | Input | Hi-Level 2 (+) |
| 5 | No Connect | |
| 6 | No Connect | |
| 7 | Input | Hi-Level 2 (-) |
| 8 | Ground | |
| 9 | Output | -DC Voltage |
| 10 | Input | Hi-Level 1 (-) |
| 11 | Ground | |
| 12 | Ground | |
| 13 | Ground | |
| 14 | Ground | |

Customer Assembly Notes

- 1) Use AWG 24 wire for all connections.
- 2) Use care to keep AWG24 wires from shorting at solder cups – keep the insulation right down to the solder cup edge.
- 3) After the wires have been soldered in place and confirmed good (using an ohm-meter), use a small amount of hot melt glue or RTV to fix the wires in place. This will help reduce stress as the cable is moved during use. Use care to not apply excess glue or the connector will be difficult to re-assemble.

| | | |
|------|-----|---|
| #1 | 1 | Connector, 14 position, gray, Redel P/N = PAGM14GLAC6SG |
| Item | Qty | Description |



TITLE: HIGH LEVEL INPUT KIT INSTRUCTION

| | | | | | |
|-----------------------------|----------------------|------------------|-------------------------------|--------------|-----------------|
| P/N(MRP #)00082180 | File # 8218A | SIZE A | DRAWING NUMBER 8218 | SCALE 1:1 | REV A |
| Dimensions are: Millimeters | Drawn By: W. Cameron | | | SHEET 1 OF 1 | |

Here's the part number: PAG.M1.4GL.AC65G

Here's the Manufacturer's webpage:

<http://www.lemousa.com/defaultwidgets.asp?nodeId=&site=454&langage=&myPage=>

Here's the webpage with the distributors:

http://intra.lemo.ch/img/distributors/LEMO_USA_Distributor_List.pdf